

PINNACLE SL300 STRAIGHT STAIRLIFT

BRINGING YOU PEACE OF MIND



DIFFERENT BY DESIGN

The Pinnacle stair lift has an established reputation as an industry leader for its innovative design, solid performance and reliability. Its patented drive system serves as a cornerstone that makes this stair lift different from all others on the market.

- Clean Operation – Pinnacle technology means no grease or lubrication. Cleaner for you, your home, or even pets.
- Efficient motor with nylon polymer gear/worm drive provides as many as up to 60 cycles on batteries during a power outage.
- Symmetrical design provides industry's most compact folded width. The rail mounts 1" from the wall and is 3" tall helping to stay out of sight.

CLEAR WALKWAYS

Harmar's exclusive automatic and mechanical folding rail (optional) removes obstructions and trip hazards at the lower landing area. Its operation is automatic, the user does not need to stop or have someone else assist. And since it's mechanical, there are no additional switches or controls.

DESIGNED WITH YOU IN MIND

Ergonomic control switch helps those with arthritis or limited hand dexterity.



Power swivel option


Seat swivels at upper landing for easy and safe entry/exit away from staircase. (power option available)

Seat belt for optimal safety

Adjustable footrest to seat height (power footrest only)

Oversized folding footrest provides more comfort and security. (power option available)

2x remote controls to call or send your stairlift

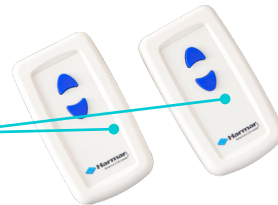

My Pinnacle gives me the confidence to go up and down my stairs safely. It operates smoothly and is very quiet. I feel comfortable in my own home and less dependent upon others. My children have peace of mind knowing I am safer with my Pinnacle.

SL300 PRODUCT SPECIFICATIONS

WEIGHT CAPACITY	300 lb	SPEED	20 FPM
TRACK LENGTH	15' 6" to 36'	INCLINE LIMITS	27° to 52°* * inclines above 45° may require local jurisdiction approval.
SEAT SIZE	19" x 16"	FOOTREST SIZE	12.6" x 14.8"
FOLDED WIDTH	- 11" from wall - 13.6" with folding rail	DRIVE	- Patented nylon polymer gear rack - Helical worm gear
COMPLIANCE	- ETL listed - CSA B44.1/ASME A17.5-2019 - Safety standard ASME A18.1-2020 - Rated residential	WARRANTY	- 2-year parts - 10-year gear rack - 1-year battery
SAFETY FEATURES	- Seat belt - Overspeed brake - Audible alert system - Swivel seat with interlock - Footrest obstruction sensor	OPTIONS	- Automatic folding rail - Chassis obstruction sensor - Extra rail + brackets - Key lock - Manual lowering tool - Power Footrest - Power Swivel - Retractable seatbelt